


Clonetics™ aortic endothelial cell systems

HAEC & D-HAEC

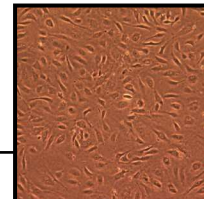

Introduction

Clonetics™ aortic endothelial cell systems contain normal human aortic endothelial cells (HAEC) or diseased human aortic endothelial cells from diabetic donors (D-HAEC) and optimized media for their growth. Each system can quickly generate HAEC or D-HAEC cultures for experimental applications in cardiovascular pharmaceutical development and vascular pathology, including atherosclerosis. Clonetics™ aortic endothelial cell systems are convenient and easy to use, allowing the researcher to focus on results. Cryopreserved HAEC and D-HAEC are shipped in third passage. Proliferating HAEC and D-HAEC are shipped in fourth passage.

Clonetics™ cells, medium and reagents are quality tested together and guaranteed to give optimum performance as a complete cell system.

Cell system components

- One aortic endothelial cell product (cryopreserved or proliferating).
- Clonetics™ EGM™-2 BulletKit™ (CC-3162) contains one 500 ml bottle of endothelial cell basal medium-2 and the following growth supplements: hydrocortisone, 0.2 ml; hFGF-B, 2 ml; VEGF, 0.5 ml; R³-IGF-1, 0.5 ml; ascorbic acid, 0.5 ml; heparin, 0.5 ml; FBS, 10 ml; hEGF, 0.5 ml; GA-1000, 0.5 ml.
- One ReagentPack™ (CC-5034) containing:

Trypsin/EDTA	100 ml
Trypsin neutralizing solution	100 ml
HEPES buffered saline solution	100 ml

Characterization of cells

Routine characterization of HAEC and D-HAEC includes immunofluorescent staining. Cells stain positive for acetylated LDL and von Willebrand (factor VIII) antigen. Cells stain negative for smooth muscle α -actin.

Performance

Recommended seeding density for subculture	2,500 - 5,000 cells/cm ²
Typical time from subculture to confluent monolayer	5 - 9 days
HAEC additional population doublings guaranteed using Clonetics™ system	15
D-HAEC additional population doublings guaranteed using Clonetics™ system	Tested through 3 passages for information only

Quality control

All cells are performance assayed and test negative for HIV-1, mycoplasma, hepatitis-B, hepatitis-C, bacteria, yeast and fungi. Cell viability, morphology and proliferative capacity are measured after recovery from cryopreservation. Clonetics™ media are formulated for optimal growth of specific types of normal human cells. Certificates of analysis (COA) for each cell strain are shipped with each order. COA for all other products are available upon request.

Lonza

Ordering information

Cryopreserved cells

CC-2535	HAEC	≥ 500,000 cells
CC-2919	D-HAEC diabetic type I	≥ 500,000 cells
CC-2920	D-HAEC diabetic type II	≥ 500,000 cells

Proliferating cells – Flasks and multiwell plates

HAEC

CC-2635		T-25 flask
CC-0222		T-75 flask
CC-0132		96-well plate

Other proliferating formats are available. Contact Scientific Support or refer to the Lonza website for details.

CC-3162	EGM™-2 BulletKit™, EBM™-2 plus SingleQuots™ of growth supplements	500 ml
CC-3156	EBM™-2, endothelial basal Medium-2	500 ml
CC-4176	EGM™-2 SingleQuots™, formulates EBM™-2 to EGM™-2	
CC-5034	ReagentPack™	
	Trypsin neutralizing solution	100 ml
	Trypsin/EDTA solution	100 ml
	HEPES buffered saline solution	100 ml

When placing an order or for technical service, please refer to the product numbers and descriptions listed above. For a complete listing of all Clonetics™ products, refer to the Lonza website or the current Lonza catalog. To obtain a catalog, additional information or technical service you may contact Lonza by web, e-mail, telephone, fax or mail.

Product Warranty

CULTURES HAVE A FINITE LIFESPAN *IN VITRO*. Lonza guarantees the performance of its cells only if Clonetics™ media and reagents are used exclusively, and the recommend protocols are followed. The performance of cells is not guaranteed if any modifications are made to the complete cell system. Cryopreserved HAEC and D-HAEC are assured to be viable and functional when thawed and maintained properly.

THESE PRODUCTS ARE FOR RESEARCH USE ONLY. Not approved for human or veterinary use, for application to humans or animals, or for use in clinical or *in vitro* procedures.

WARNING: CLONETICS™ AND POIETICS™ PRODUCTS CONTAIN HUMAN SOURCE MATERIAL, TREAT AS POTENTIALLY INFECTIOUS. Each donor is tested and found non-reactive by an FDA approved method for the presence of HIV-I, hepatitis B Virus and hepatitis C virus. Where donor testing is not possible, cell products are tested for the presence of viral nucleic acid from HIV, hepatitis B Virus, and hepatitis C virus. Testing can not offer complete assurance that HIV-1, hepatitis B virus, and hepatitis C virus are absent. All human sourced products should be handled at the biological safety level 2 to minimize exposure of potentially infectious products, as recommended in the CDC-NIH manual, [Biosafety in Microbiological and Biomedical Laboratories](#), 5th edition. If you require further information, please contact your site safety officer or Scientific Support.